

HYPROP



Generate Moisture Release Curves


Forget the pain of using pressure plates. Use Hyprop to generate detailed moisture release curves, typically within a day.

Just take your undisturbed sample, insert the measurement head, put it on the balance, and one to two days, you have a detailed moisture release curve.

- Uses undisturbed soil samples.*
- Works in all soil types*
- Bonus feature gives unsaturated hydraulic conductivity values for the soil sample.*
- Range: 0 to -85 kPa*

Hyprop Generates Moisture Release Curves Automatically

Drying portion of the moisture curve generated by Hyprop.

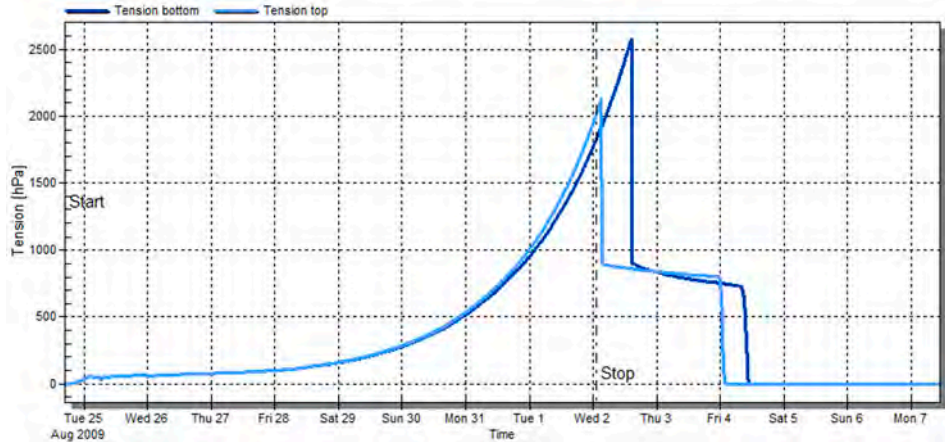

Hyprop Specifications

Range	+20 hPa to -1200 hPa / -2500 hPa
Resolution	0.01 hPa
Accuracy	1.5 hPa (0hPa to 800 hPa)
Interface	RS485 tensioLINK

Hyprop w/ Laboratory scale (optional)

Range	0 to 2.5kg
Resolution	0.01g
Accuracy	0.1g
Interface	RS232



Optional Refilling kit, makes using the Hyprop easier.

© 2009 DECAGON
Printed in USA



2365 NE Hopkins Court
Pullman, Washington 99163
1-800-755-2751
www.decagon.com